


OMRON®

INSTRUCTION MANUAL

5 series Blood Pressure Monitor Model BP742

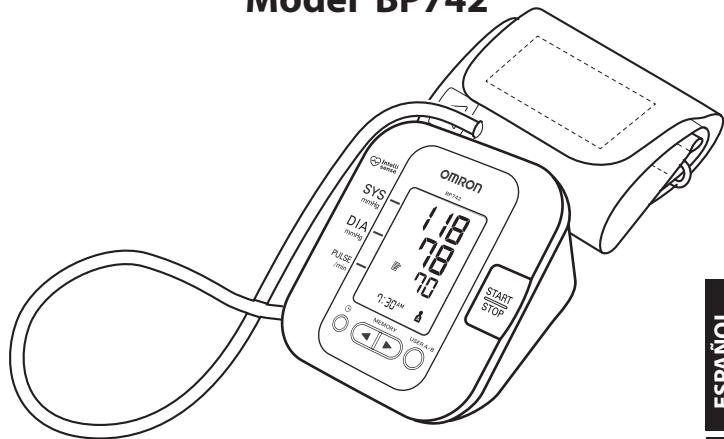

ESPAÑOL

ENGLISH

TABLE OF CONTENTS

Before using the Monitor

Introduction	3
Safety Information	4
Operating the Device	4
Care and Maintenance	5

Operating Instructions

Before Taking a Measurement	6
Know Your Unit	7
Unit Display	7
Display Symbols	8
Irregular Heartbeat Symbol	8
Movement Error Symbol	8
Heartbeat Symbol/Hypertension Indicator	8
Battery Installation	10
Setting the Date and Time	12
Applying the Arm Cuff	14
Taking a Measurement	16
Special Conditions	18
Using the Memory Function	19
To Display the Measurement Values	19
To Delete All Values Stored in the Memory	20

Care and Maintenance	21
Error Indicators	22
Troubleshooting Tips	23
FCC Statement	24
Limited Warranties	25
Specifications	26

INTRODUCTION

Thank you for purchasing the OMRON® BP742 IntelliSense® Automatic Blood Pressure Monitor.

Fill in for future reference.

DATE PURCHASED: _____

SERIAL NUMBER: _____

- *Staple your purchase receipt here*
- *Register your product on-line at www.register-omron.com*

Your new digital blood pressure monitor uses the oscillometric method of blood pressure measurement. This means the monitor detects your blood's movement through your brachial artery and converts the movements into a digital reading. An oscillometric monitor does not need a stethoscope so the monitor is simple to use.

The BP742 comes with the following components:

- Monitor • Standard Adult Arm Cuff • 4 "AA" Alkaline Batteries
- Instruction Manual printed in English and Spanish • Quick Start Guide

The BP742 IntelliSense® Automatic Blood Pressure Monitor is intended for home use.

SAVE THESE INSTRUCTIONS

SAFETY INFORMATION

To assure the correct use of the product, basic safety measures should always be followed including the warnings and cautions listed in this instruction manual.

SAFETY SYMBOLS USED IN THIS INSTRUCTION MANUAL	
⚠ WARNING	Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.
⚠ CAUTION	Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user or patient or damage to the equipment or other property.

OPERATING THE DEVICE

- ⚠ Contact your physician for specific information about your blood pressure. Self-diagnosis and treatment using measured results may be dangerous. Follow the instructions of your physician or licensed healthcare provider.
- ⚠ The monitor is not intended to be a diagnostic device.
- ⚠ If battery fluid should get in your eyes, immediately rinse with plenty of clean water. Contact a physician immediately.
- ⚠ Read all of the information in the instruction manual and any other literature in the box before operating the unit.
- ⚠ This device is intended for use in measuring blood pressure and pulse rate in the adult population. Do not use this device on infants or persons who cannot express their intentions.
- ⚠ Read the Special Conditions section of the instruction manual if your systolic pressure is known to be more than 220 mmHg. Inflating to a higher pressure than necessary may result in bruising where the cuff is applied.
- ⚠ Operate the device only as intended. Do not use the device for any other purpose.
- ⚠ Dispose of the device, components and optional accessories according to applicable local regulations. Unlawful disposal may cause environmental pollution.

SAFETY INFORMATION

- ⚠ Do not use a cellular phone near the device. It may result in an operational failure.
- ⚠ Use only Omron authorized parts and accessories. Parts and accessories not approved for use with the device may damage the unit.
- ⚠ Use only 1.5 V “AA” alkaline batteries with this device. Do not use other types of batteries. This may damage the unit

CARE AND MAINTENANCE

- ⚠ Do not subject the monitor to strong shocks, such as dropping the unit on the floor.
- ⚠ Do not submerge the device or any of the components in water.
- ⚠ Store the device and the components in a clean, safe location.
- ⚠ Changes or modification not approved by the manufacturer will void the user warranty. Do not disassemble or attempt to repair the unit or components.


BEFORE TAKING A MEASUREMENT

To ensure a reliable reading follow these recommendations:


1. Avoid eating, drinking alcohol, smoking, exercising, and bathing for 30 minutes before taking a measurement. Rest for at least 15 minutes before taking the measurement.
2. Stress raises blood pressure. Avoid taking measurements during stressful times.
3. Measurements should be taken in a quiet place.
4. Remove tight-fitting clothing from your left arm.
5. Sit on a chair with your feet flat on the floor. Rest your left arm on table so that the cuff is at the same level as your heart.
6. Remain still and do not talk during the measurement.
7. Wait 2-3 minutes between measurements. The wait time allows the arteries to return to the condition prior to taking the blood pressure measurement. You may need to increase the wait time depending on your individual physiological characteristics.
8. Keep a record of your blood pressure and pulse readings for your physician. A single measurement does not provide an accurate indication of your true blood pressure. You need to take and record several readings over a period of time. Try to measure your blood pressure at the same time each day for consistency.

KNOW YOUR UNIT

Main Unit:

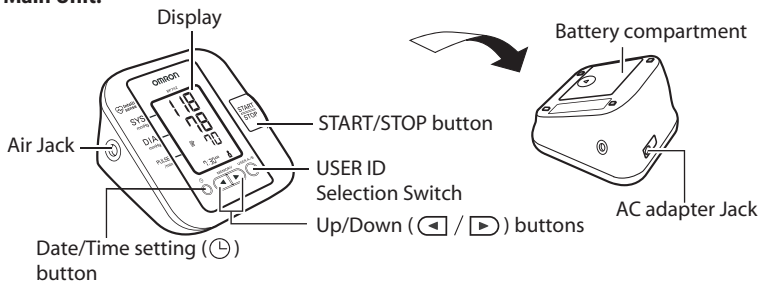

Components:

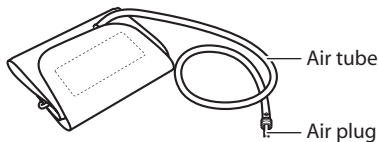

Standard Adult Arm Cuff




Instruction Manual




Quick Start Guide

UNIT DISPLAY

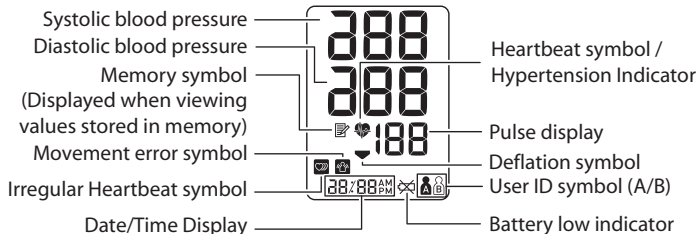




DISPLAY SYMBOLS

IRREGULAR HEARTBEAT SYMBOL ()


When the monitor detects an irregular rhythm two or more times during the measurement, the Irregular Heartbeat Symbol () will appear on the display with the measurement values.

An irregular heartbeat rhythm is defined as a rhythm that is more than 25% slower or 25% faster from the average rhythm detected while the monitor is measuring the systolic and diastolic blood pressure.

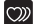If the irregular heartbeat symbol () displays with your measurement results, we recommend you consult your physician. Follow the directions of your doctor.

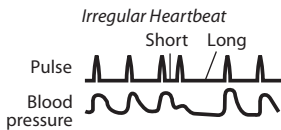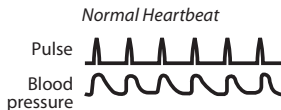




MOVEMENT ERROR SYMBOL ()

The Movement Error Symbol () is displayed if you move your body during the measurement. Please remove the arm cuff, and wait 2-3 minutes. Reapply the arm cuff and take another measurement.


HEARTBEAT SYMBOL/HYPERTENSION INDICATOR ()

TAKING A MEASUREMENT

The Heartbeat Symbol () flashes on the display during the measurement. The Heartbeat Symbol () flashes at every heartbeat.

When the measurement is complete, the Heartbeat Symbol () flashes on the display with your blood pressure and pulse rate if the reading is above 135 for the Systolic Blood Pressure value and/or above 85 for the Diastolic Blood Pressure value.

USING THE MEMORY FUNCTION

The Heartbeat Symbol () flashes when the blood pressure reading stored in the memory is above 135 for the Systolic Blood Pressure value and/or above 85 for the Diastolic Blood Pressure value.

DISPLAY SYMBOLS

The American Heart Association recommends the following guideline for upper limit of normal for home blood pressure readings:

AHA Home Guideline for Upper Limit of Normal BP	
Systolic Blood Pressure	135 mmHg
Diastolic Blood Pressure	85 mmHg


This is a general guideline because blood pressure depends on age and morbidity. The AHA recommends a lower home blood pressure goal for certain patients, including diabetic patients, pregnant women, and patients with renal failure¹.

⚠ WARNING

The monitor is not intended to be a diagnostic device.

⚠ WARNING

Contact your physician for specific information about your blood pressure. Self-diagnosis and treatment using measured results may be dangerous. Follow the instructions of your physician or licensed healthcare provider.

Omron Healthcare recommends you contact your physician if the Heartbeat Symbol () flashes on the display with your measurement values.



¹American Heart Association Scientific Statement 2004

BATTERY INSTALLATION

⚠ CAUTION

Use only 1.5V “AA” alkaline batteries with this device.

1. Press the ▲ indicator on the battery cover and slide the cover off in the direction of the arrow.
2. Install 4 “AA” size batteries so the + (positive) and - (negative) polarities match the polarities of the battery compartment as indicated.
3. Replace the battery cover.

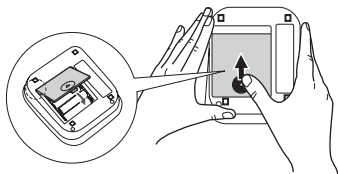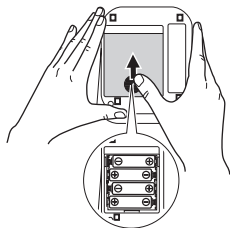

BATTERY INSTALLATION

BATTERY REPLACEMENT



Low Battery Indicator

When the Low Battery Indicator appears on the display, turn the monitor off and remove all the batteries. Replace with four new batteries at the same time. Long-life alkaline batteries are recommended.

⚠ WARNING

If battery fluid should get in your eyes, immediately rinse with plenty of clean water. Contact a physician immediately.

- NOTES:**
- Turn the unit off before replacing the batteries. If the batteries are removed while the unit is still on, the date and time will be reset to that of the previous use. However, the measurement results are not deleted.
 - When the batteries are replaced, you may need to reset the date and time. If “-:-” appears on the display, refer to “Setting the Date and Time”.
 - The supplied batteries may have a shorter life.

⚠ CAUTION


Dispose of the device, components and optional accessories according to applicable local regulations. Unlawful disposal may cause environmental pollution.

SETTING THE DATE AND TIME


Set the monitor to the current date and time before taking a measurement for the first time. If the date and time has not been set, the measurement values are stored without any date or time information.

1. TO START THE SETTING

Press the Date/Time Setting button (⊖) to adjust the date and time.

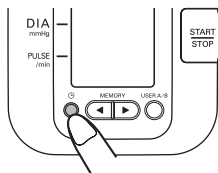

The year flashes on the display

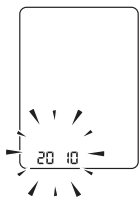


2. SETTING THE YEAR


The year can be set between 2010 and 2040. When the display reaches 2040, it will return to 2010.

Press the Up button (▶) to advance by increments of one year.

Press the Date/Time Setting button (⊖) to set the current year.

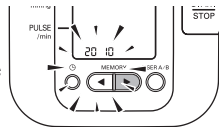

The month flashes on the display.

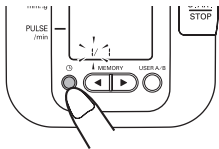

- NOTES:**
- Press the Down button (◀) to decrease the date and time values by increments of one.
 - Press and hold the button to increase or decrease the date and time values faster.


SETTING THE DATE AND TIME

3. SETTING THE MONTH

Press the Up button (▶) to advance by increments of one month.

Press the Date/Time Setting button (⌚) to set the current month.

The day flashes on the display.

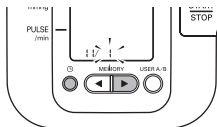


4. SETTING THE DAY

Press the Up button (▶) to advance by increments of one day.

Press the Date/Time Setting button (⌚) to set the current day.

The hour flashes on the display.

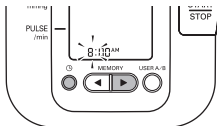


5. SETTING THE HOUR

The time is set using AM or PM.

Press the Up button (▶) to advance by increments of one hour.

Press the Date/Time Setting button (⌚) to set the current hour.

The minute flashes on the display.

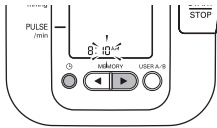

6. SETTING THE MINUTE

Press the Up button (▶) to advance by increments of one minute.

Press the Date/Time Setting button (⌚) to set the minute setting.

7. TO TURN THE MONITOR OFF

Press the START/STOP button.

NOTE: The monitor will automatically turn off five minutes after the minute has been set.

APPLYING THE ARM CUFF

⚠ CAUTION


Read all of the information in the instruction manual and any other literature in the box before operating the unit.

⚠ CAUTION

This device is intended for use in measuring blood pressure and pulse rate in the adult population. Do not use this device on infants or persons who cannot express their intentions.


APPLYING THE CUFF ON THE LEFT ARM

1. Make sure the air plug is securely inserted in the main unit.
2. Remove tight-fitting clothing from your upper arm.

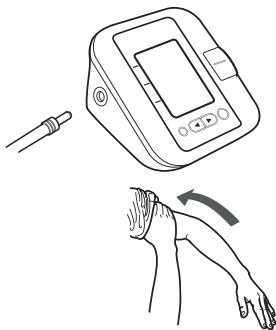

APPLYING THE ARM CUFF


3. Apply the cuff to your left upper arm so the arrow is centered on the inside of your arm and aligned with your middle finger. The air tube runs down the inside of your arm. The bottom of the cuff should be approximately 1/2 " above your elbow.
4. Wrap the cuff firmly in place around your arm using the cloth fastening.

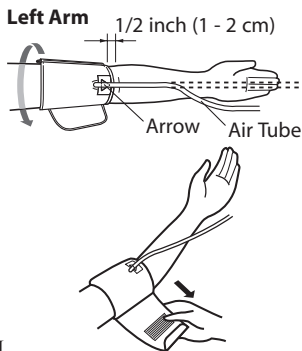

APPLYING THE CUFF ON THE RIGHT ARM

When taking a measurement using the right arm, use this instruction for Step 3.

Apply the cuff to your right upper arm so the arrow and the air tube runs along the inside of your arm. The bottom of the cuff should be approximately 1/2 " above your elbow.

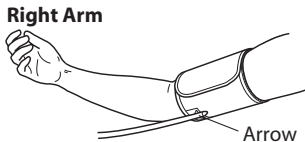


NOTE: Be careful not to rest your arm on the air tube. This will restrict the flow of air to the cuff.


TAKING A MEASUREMENT

1. Press the START/STOP button to turn the monitor on.

All symbols appear on the display.

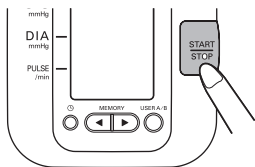

2. Select your USER ID (A or B).


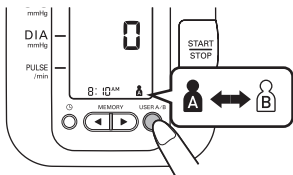


3. Sit on a chair with your feet flat on the floor.
Place your arm on a table so the cuff is level with your heart.
Keep still and do not talk during measurement.
4. Press the START/STOP button.



The cuff starts to inflate automatically. As the cuff inflates, the monitor automatically determines your ideal inflation level.


Remain still and do not move your arm until the entire measurement process is completed.

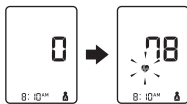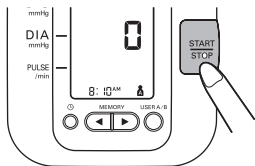

NOTE: To stop the inflation or measurement, push the START/STOP button. The monitor will stop inflating, start deflating, and will turn off.

TAKING A MEASUREMENT


5. Inflation stops automatically and the measurement is started. As the cuff deflates, decreasing numbers appear on the display and the Heartbeat Symbol (♥) flashes.

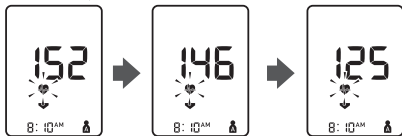

6. When the measurement is complete, the arm cuff completely deflates. Your blood pressure and pulse rate are displayed.



7. Press the START/STOP button to turn the monitor off.

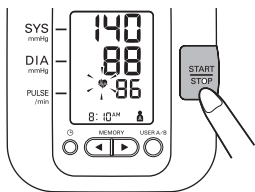

NOTE: The monitor will automatically turn off after five minutes.


TAKING A MEASUREMENT

SPECIAL CONDITIONS

⚠ CAUTION

Inflating to a higher pressure than necessary may result in bruising where the cuff is applied.

If your systolic pressure is known to be more than 220 mmHg, press and hold the START/STOP button until the monitor inflates 30 to 40 mmHg higher than your expected systolic pressure.

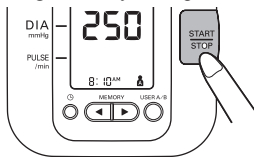

- NOTES:**
- Do not apply more pressure than necessary.
 - The monitor will not inflate above 299 mmHg.

When the measurement is complete, the arm cuff completely deflates. Your blood pressure and pulse rate are displayed.

USING THE MEMORY FUNCTION


The monitor is designed to store the blood pressure and the pulse rate in the memory for two people (USER A and USER B) every time a measurement is completed. The monitor automatically stores up to 30 sets of measurement values (blood pressure and pulse rate) for each user (A and B). When 30 sets of measurement values are stored, the oldest record is deleted to save the most recent values.

TO DISPLAY THE MEASUREMENT VALUES

1. Press the START/STOP button to turn the monitor on.
2. Select your USER ID (A or B).
3. Press the Up button (▶) to display the most recent set of measurement values on the screen.
Press the Down button (◀) to display the oldest set of measurement values.

USING THE MEMORY FUNCTION

- NOTES:**
- The measurement value number is displayed briefly in the Pulse Display area when you press the Up/Down buttons (◀ / ▶).
 - The date and time are alternately displayed with the measurement values.
 - If the date and time are not set, “-:--” is displayed with the measurement values.
 - Press the button repeatedly to display the next values.
 - Press and hold the button to display the values faster.

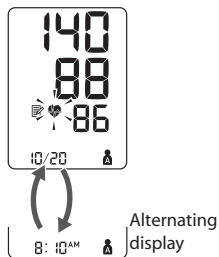


TO DELETE ALL VALUES STORED IN THE MEMORY

The values stored in the memory are deleted by USER ID.

You cannot partially delete values stored in the memory. All values for the user you select will be deleted.

1. Press the START/STOP button to turn the monitor on.
2. Select your USER ID (A or B).
3. While holding the Down button (◀), press the START/STOP button down for more than five seconds until the display appears as it does in the following illustration.

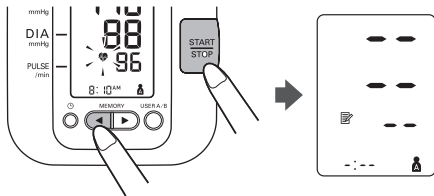

4. Press the START/STOP button to turn the monitor off.

NOTE: The monitor will automatically turn off after five minutes.

CARE AND MAINTENANCE

To keep your digital blood pressure monitor in the best condition and protect the unit from damage, follow the directions listed below:

Do not forcefully bend the arm cuff or air tube. Do not fold tightly.

Clean the monitor with a soft dry cloth.

Do not use any abrasive or volatile cleaners.

Do not attempt to clean the cuff.

⚠ CAUTION

Do not submerge the device or any of the components in water. Do not subject the monitor to extreme hot or cold temperatures, humidity or direct sunlight.

⚠ CAUTION

Store the device and the components in a clean, safe location.

⚠ CAUTION

Do not subject the monitor to strong shocks, such as dropping the unit on the floor.

Remove the batteries if the unit will not be used for three months or longer. Always replace all the batteries with new ones at the same time.

Use the unit consistent with the instruction provided in this manual.

Use only authorized parts and accessories.

Use the unit in a manner consistent with the instructions provided in this manual.








⚠ CAUTION

Use only Omron authorized parts and accessories. Parts and accessories not approved for use with the device may damage the unit.

⚠ CAUTION

Changes or modification not approved by the manufacturer will void the user warranty. Do not disassemble or attempt to repair the unit or components.

ERROR INDICATORS

SYMBOL	CAUSE	CORRECTION
	Monitor could not detect pulse wave.	Take another measurement and remain still until the measurement is complete.
	Cuff under-inflated. Cuff not applied correctly.	Remove the arm cuff. Read “Taking a Measurement”. Wait 2-3 minutes. Take another measurement.
	Cuff over-inflated.	
	Device error.	Contact Customer Service.
	Movement during measurement.	Keep still and take another measurement.
	Irregular or weak pulses are detected.	Remove the arm cuff. Wait 2-3 minutes and then take another measurement. Repeat the steps in section “Applying the Arm Cuff.” If this error continues to appear, contact your doctor.
	Batteries are exhausted.	Replace the four batteries. Refer to “Battery Installation.”

TROUBLESHOOTING TIPS

PROBLEM	CAUSES AND SOLUTIONS
No power. No display appears on the unit.	Replace all four batteries with new ones. Check the battery installation for proper placement of the battery polarities.
Measurement values appear too high or too low.	Blood pressure varies constantly. Many factors including stress, time of day, how you wrap the cuff, may affect your blood pressure. Review the sections “Before Taking a Measurement” and “Taking a Measurement”.

FCC STATEMENT

Note:

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE (for U.S.A. only)

This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. The product generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the product on and off, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and the receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

POTENTIAL FOR RADIO/TELEVISION INTERFERENCE (for Canada only)

This digital apparatus does not exceed the Class B limits for radio noise emissions from digital apparatus as set out in the interference-causing equipment standard entitled “Digital Apparatus”, ICES-003 of the Canadian Department of Communications.

Cet appareil numérique respecte les limites de bruits radioélectriques applicables aux appareils numériques de Classe B prescrites dans la norme sur le matériel brouilleur: “Appareils Numériques”, ICES-003 édictée par le ministre des communications.

Changes or modifications not expressly approved by the party responsible for compliance could void the user’s authority to operate the equipment.

LIMITED WARRANTIES

Your BP742 IntelliSense® Automatic Blood Pressure Monitor, excluding the arm cuff, is warranted to be free from defects in materials and workmanship appearing within 2 years from the date of purchase, when used in accordance with the instructions provided with the monitor. The arm cuff is warranted to be free from defects in materials and workmanship appearing within one year from the date of purchase when the monitor is used in accordance with the instructions provided with the monitor. The above warranties extend only to the original retail purchaser.

We will, at our option, repair or replace without charge any monitor or arm cuff covered by the above warranties. Repair or replacement is our only responsibility and your only remedy under the above warranties.

To obtain warranty service contact Customer Service by calling **1-800-634-4350** for the address of the repair location and the return shipping and handling fee.

Enclose the Proof of Purchase. Include a letter, with your name, address, phone number, and description of the specific problem. Pack the product carefully to prevent damage in transit. Because of possible loss in transit, we recommend insuring the product with return receipt requested.

THE FOREGOING IS THE SOLE WARRANTY PROVIDED BY OMRON IN CONNECTION WITH THIS PRODUCT, AND OMRON HEREBY DISCLAIMS ANY OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. IMPLIED WARRANTIES AND OTHER TERMS THAT MAY BE IMPOSED BY LAW, IF ANY, ARE LIMITED IN DURATION TO THE PERIOD OF THE ABOVE EXPRESS WARRANTY.

OMRON SHALL NOT BE LIABLE FOR LOSS OF USE OR ANY OTHER SPECIAL, INCIDENTAL, CONSEQUENTIAL OR INDIRECT COSTS, EXPENSES OR DAMAGES.

This warranty provides you with specific legal rights, and you may have other rights that vary by jurisdiction. Because of special local requirements, some of the above limitations and exclusions may not apply to you.

FOR CUSTOMER SERVICE

Visit our web site at:
Call toll free:

www.omronhealthcare.com
1-800-634-4350

SPECIFICATIONS

Model:	BP742 (HEM-7200-Z)
Display:	LCD Digital Display
Measurement Range:	Pressure: 0 to 299 mmHg, Pulse: 40 to 180/min.
Accuracy/Calibration:	Pressure: ± 3 mmHg or 2% of reading Pulse: $\pm 5\%$ of reading
Inflation:	Automatic by electric pump
Deflation:	Automatic pressure release valve
Rapid Air Release:	Automatic exhaust valve
Measurement Method:	Oscillometric method
Power Source:	1.5V 4 "AA" batteries
Battery Life:	Approx. 300 uses with 4 new alkaline batteries
Operating Temperature /Humidity:	50°F to 104°F (10°C to 40°C) /15 to 90% RH
Storage Temperature /Humidity/Air Pressure:	-4°F to 140°F (-20°C to 60°C) /10 to 95% RH 700-1060 hPa
Main Unit Weight:	Approximately 12 oz. (340g) not including batteries
Main Unit Dimensions:	Approx. 5 1/2"(l) \times 4 7/8"(w) \times 3 3/8" (h) (141 mm \times 123 mm \times 85 mm)
Cuff Size:	Approximately 5 3/4"(w) \times 17 1/2"(l) (146 mm \times 446 mm) Cuff tube 23 5/8" (600 mm)
Cuff Circumference:	Fits arm circumferences 9" to 13" (220 mm to 320 mm)
Memory:	Up to 30 readings per user
Contents:	Main Unit, Standard Adult Arm Cuff, 4 "AA" Alkaline Batteries, Instruction Manual, and Quick Start Guide

SPECIFICATIONS

Optional Accessories: Large Cuff 13" - 17" (320 mm - 420 mm), AC Adapter



= Type B

NOTE: These specifications are subject to change without notice.

Omron Healthcare, Inc.
1200 Lakeside Drive
Bannockburn, IL 60015
www.omronhealthcare.com
Copyright © 2010 Omron Healthcare, Inc.

Made in China

5329417-5A